

Autotome™ RX

Cannulating Sphincterotome



RX LOCKING DEVICE†

Guidewire locks in place to maintain access

- ◆ Accepts short 260cm guidewire or a standard 450cm guidewire to meet physician technique preferences
- ◆ Guidewire may be locked and unlocked at any time during the procedure
- ◆ Facilitates device exchange, contrast injection and biliary therapy
- ◆ RX 39 compatible with .025" guidewire, RX 44 and RX 49 compatible with .035" guidewire



AUTOTOME RX CANNULATING SPHINCTEROTOME*

Product Number	Description	Cut Wire Length (mm)	Tip Length (mm)	Tip O.D. (Fr)	Recommended Guidewire
M00545150	Autotome RX 49	20	5	4.9	.035" Hydra Jagwire™ Guidewire
M00545160	Autotome RX 49	30	5	4.9	.035" Hydra Jagwire Guidewire
M00545170	Autotome RX 44	20	5	4.4	.035" Hydra Jagwire Guidewire
M00545180	Autotome RX 44	30	5	4.4	.035" Hydra Jagwire Guidewire
M00545190	Autotome RX 39	20	5	3.9	.025" Jagwire™ Guidewire
M00545200	Autotome RX 39	30	5	3.9	.025" Jagwire Guidewire

*Compatible with Universal Active Cords.

European Headquarters – Paris

T: +33 1 57 66 80 00 F: +33 1 57 66 84 99

Argentina T: +5411 4896 8556 F: +5411 4896 8550

Australia/New Zealand T: +61 2 8336 5555 F: +61 2 8336 0404

Austria T: +43 1 726 30 05 1022 F: +43 1 726 30 05 1000

Balkans T: +30 210 95 37 890 F: +30 210 95 79 836

Belgium (Freephone) T: 080094 494 F: 080093 343

Brazil T: +55 11 5502 8500 F: +55 11 5103 2212

Bulgaria T: +359 2 986 50 48 F: +359 2 986 57 09

Canada T: +1 888 359 9691 F: +1 888 575 7396

Chile T: +562 445 4904 F: +562 445 4915

China – Beijing T: +86 10 8525 1588 F: +86 10 8525 1566

China – Guangzhou T: +86 20 8767 9791 F: +86 20 8767 9789

China – Shanghai

T: +86 21 6141 5959 F: +86 21 6141 5900

Colombia T: +57 1 629 5045 F: +57 1 612 4761

Czech Republic T: +420 2 3536 2911 F: +420 2 3536 4334

Denmark (Freephone) T: 80 30 80 05 F: 80 30 80 05

Finland (Freephone) T: 09 622 42 85 F: 09 622 42 86

France T: +33 1 39 30 49 30 F: +33 1 39 30 49 31

Germany (Freephone) T: 0800 072 3303 F: 0800 072 3319

Greece T: +30 210 95 37 890 F: +30 210 95 79 836

Hong Kong T: +852 2960 7100 F: +852 2563 5276

Hungary T: +36 1 456 30 40 F: +36 1 456 30 41

India – Bangalore T: +91 80 5112 11045 F: +91 80 5112 1106

India – Chennai T: +91 44 2648 0318 F: +91 44 2641 4695

India – Delhi

T: +91 11 2618 0445/6 F: +91 11 2618 1024

India – Mumbai T: +91 22 5677 8844 F: +91 22 2617 2783

Italy T: +39 010 60 60 1 F: +39 010 60 60 200

Korea T: +82 2 3476 2121 F: +82 2 3476 1776

Mexico T: +52 55 5687 63 90 F: +52 55 5687 62 28

Middle East / Gulf / North Africa T: +961 1 805 282 F: +961 1 805 445

The Netherlands T: +31 43 356 5260 F: +31 43 356 6265

Norway (Freephone) T: 800 104 04 F: 800 101 90

Poland T: +48 22 435 1414 F: +48 22 435 1410

Portugal T: +351 1 381 25 40 F: +351 1 381 28 29

South Africa T: +27 11 840 8800 F: +27 11 463 6077

South East Asia – Malaysia T: +60 3 2284 1459 F: +60 3 2284 1507

South East Asia – Philippines

T: +63 2 687 3239 F: +63 2 687 3047

South East Asia – Singapore T: +65 6418 8888 F: +65 6418 8889

South East Asia – Thailand T: +66 2 2654 3810 F: +66 2 2654 3818

Spain T: +34 901 11 12 15 F: +34 93 405 9045

Sweden (Freephone) T: 020 65 25 30 F: 020 65 25 35

Switzerland (Freephone) T: 0800 826 796 F: 0800 826 787

Taiwan T: +886 2 2747 7278 F: +886 2 2747 7270

Turkey – Istanbul T: +90 216 464 36 66 F: +90 216 464 36 67

Uruguay T: +59 82 900 6212 F: +59 82 900 6212

UK & Eire T: +44 1227 866633 F: +44 1227 830371

Venezuela T: +58 212 959 8106 F: +58 212 959 828

Boston Scientific

Delivering what's next.™

www.bostonscientific.com/international
www.bostonscientific.ca
www.bostonscientific.co.uk
www.bostonscientific.de
www.bostonscientific.es
www.bostonscientific.fr
www.bostonscientific.ie
www.bostonscientific.it
www.bostonscientific.nl

© 2006 Boston Scientific Corporation or its affiliates. All rights reserved.
DINEND2152EA

Autotome™ RX

Cannulating Sphincterotome



Boston Scientific

RX Biliary System™
A difference you can feel.™

†RX Locking Device sold separately.

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. Information for the use only in countries with applicable health authority product registrations.
PSST IC 1827 Printed in the UK by Gosling.

RX Biliary System™

A difference you can feel.™

The endoscope shows you the way, but it's your sense of touch that gets you where you need to go. The devices that make up the RX Biliary System have been developed to provide access and control with a sensation of device responsiveness that truly is... a difference you can feel.

Over the years, Boston Scientific has earned the confidence and allegiance of customers around the world who have adopted the RX Biliary System for everyday use. Our legacy of success in developing high performance RX biliary devices has enabled us to change the landscape of ERCP procedures. By innovating and refining products, we've put control and access right at your fingertips.

Autotome™ RX Cannulating Sphincterotome

FLEXIBLE AND RESPONSIVE DESIGN

Direct Wire Technology™ is designed for flexibility and responsiveness to aid in precise sphincterotomy

- ◆ Tip orientation and cut-wire positioning capabilities developed to facilitate cannulation
- ◆ Range of tapered tip sizes: 3.9Fr, 4.4Fr and 4.9Fr tip sizes to accommodate varying patient anatomy
- ◆ 5cm – 7cm exchange length engineered for single-operator control
- ◆ Option to transition from over-the-wire to short-wire ERCP technique without removing devices from the endoscope

MERGING LUMEN TECHNOLOGY™

Proximal lumens merge into a single distal lumen on RX 39 and RX 44 Sphincterotomes

- ◆ Developed to facilitate contrast injection whilst maintaining guidewire access through open lumen technology
- ◆ May help reduce the risk of submucosal injection by allowing contrast flow to be diverted



RX 39



RX 44




RX 49



TIP ORIENTATION

Ergonomic, rotating handle engineered for excellent control

- ◆ Designed to aid precise tip movement and wire manipulation to facilitate cannulation
- ◆ Handles labelled with respective tip outer diameter for ease of selection*

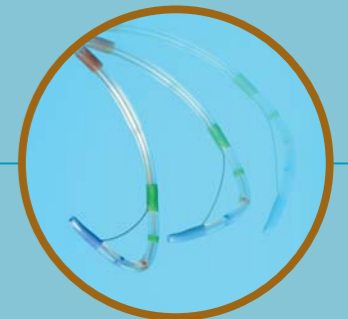

FLEXCONTROL SYSTEM*

Developed to aid in precise sphincterotomy

- ◆ Engineered to enhance guidewire manipulation and bowing control

*Autotome RX 39 and RX 44 Sphincterotomes only.