


Interject™

Sclerotherapy Needles

Boston
Scientific

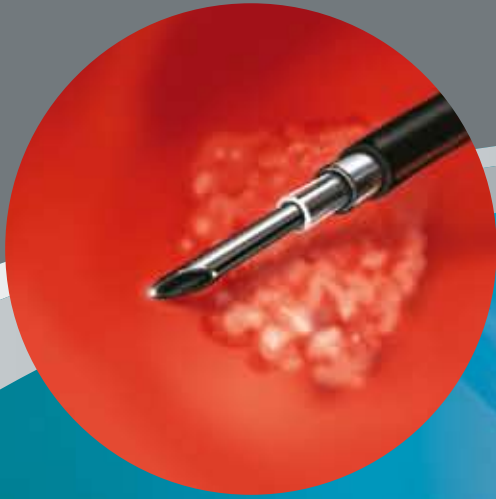

Performance
You can count on



Interject™

Sclerotherapy Needles

Designed to provide **performance you can count on** for injection therapy applications.

TORTUOUS ANATOMY

Unique star-shaped inner catheter design* helps promote reliable needle actuation in tortuous anatomy (even after kinking).

EFFICIENT ACCESS

Strong, pushable sheath construction designed to facilitate advancement through difficult pathways to remote anatomy.

SMOOTH NEEDLE PENETRATION


Beveled needle designed to facilitate injection.

SECURE NEEDLE POSITIONING

Red Spacer Clip is designed to maintain needle position within the sheath during scope passage.

RANGE OF OPTIONS

- 200cm and 240cm sheath lengths
- 4mm and 6mm needle lengths
- 23 and 25 gauge needles
- Clear and contrast sheath available.

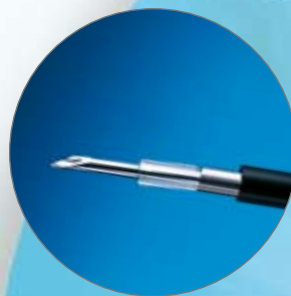

Smooth beveled needle designed to facilitate injection



Packaged with new Red Spacer Clip which is designed to keep the needle inside the catheter until used



Convenient locking hub for procedural efficiency

ORDERING INFORMATION

Order Number	Sheath O.D. (mm)	Sheath Length (cm)	Sheath Design	Needle Gauge	Needle O.D. (mm)	Needle I.D. (mm)	Maximum Needle Extension Length (mm)
M00518151	1.8	200	Contrast	23	0.64	0.32	4
M00518161	1.8	200	Contrast	25	0.51	0.24	4
M00518251	1.8	200	Contrast	23	0.64	0.32	6
M00518261	1.8	200	Contrast	25	0.51	0.24	6
M00518351	2.3	240	Contrast	23	0.64	0.32	4
M00518361	2.3	240	Contrast	25	0.51	0.24	4
M00518111	1.8	200	Clear	25	0.51	0.24	4
M00518301	2.3	240	Clear	23	0.64	0.32	4
M00518311	2.3	240	Clear	25	0.51	0.24	4

Packaged 5 Per Box

*Patent on file Boston Scientific Corp.

European Distribution Centre – The Netherlands
T: +31 45 54 67 700 F: +31 45 54 67 800

European Headquarters – Paris
T: +33 1 57 66 80 00 F: +33 1 57 66 84 99

Argentina (Freephone)
T: 0800 555 2678 F: +5411 4896 8517

Australia/New Zealand
T: +61 2 8063 8100 F: +61 2 9330 1404

Austria
T: +43 1 60 81 810 F: +43 1 60 810 60

Belgium (Freephone)
T: 0800 94 494 F: 0800 93 343

Brazil
T: +55 11 5502 8500 F: +55 11 5103 2212

Canada
T: +1 888 359 9691 F: +1 888 575 7396

Chile
T: +562 445 4904 F: +562 445 4915

China – Beijing
T: +86 10 8525 1588 F: +86 10 8525 1566

China – Guangzhou
T: +86 20 8767 9791 F: +86 20 8767 9789

China – Shanghai
T: +86 21 6141 5959 F: +86 21 6141 5900

Colombia
T: +57 1 629 5045 F: +57 1 612 4761

Czech Republic
T: +420 296 331 901 F: +420 296 331 935

Denmark (Freephone)
T: 80 30 80 02 F: 80 30 80 05

Eire (Freephone)
T: 1800 882 969 F: 1800 882 968

Finland
T: +358 20 762 8882 F: +358 20 762 8883

France
T: +33 1 39 30 49 00 F: +33 1 39 30 49 01

Germany (Freephone)
T: 08000 723300 F: 08000 723319

Greece
T: +30 210 9542 300 F: +30 210 9542 310

Hong Kong
T: +852 2960 7100 F: +852 2563 5276

Hungary
T: +36 1 456 30 40 F: +36 1 456 30 41

India – Bangalore
T: +91 80 2212 4928/9 F: +91 80 2207 5153

India – Chennai
T: +91 44 2220 1879 F: +91 44 2220 0716

India – Delhi
T: +91 11 4243 2222 F: +91 11 2610 0808

India – Mumbai
T: +91 22 4030 9165 F: +91 22 4040 9199

Italy
T: +39 010 60 60 1 F: +39 010 60 60 200

Korea
T: +82 2 3476 2121 F: +82 2 3476 1776

Mexico
T: +52 55 5687 63 90 F: +52 55 5687 62 28

Middle East/Gulf/North Africa
T: +961 1 805 410 F: +961 1 805 445

The Netherlands
T: +31 30 602 55 44 F: +31 30 602 55 05

Norway (Freephone)
T: 800 104 04 F: 800 101 90

Poland
T: +48 22 435 14 14 F: +48 22 435 14 10

Portugal
T: +351 1 381 25 40 F: +351 21 381 25 58

South Africa
T: +27 11 840 8600 F: +27 11 463 6077

South East Asia – Malaysia
T: +60 3 2283 3813 F: +60 3 2284 3813

South East Asia – Philippines
T: +63 2 687 6994 F: +63 2 687 3047

South East Asia – Singapore
T: +65 6418 8888 F: +65 6418 8899

South East Asia – Thailand
T: +66 2 6380 100 F: +66 2 6380 400

South East Europe
T: +30 210 9542 300 F: +30 210 9542 310

Spain – Barcelona
T: +34 93 444 72 00 F: +34 93 405 90 45

Spain – Madrid
T: +34 901 11 12 15 F: +34 91 319 50 03

Sweden
T: +46 42 25 69 00 F: +46 42 25 69 69

Switzerland
T: +41 32 626 57 00 F: +41 32 626 57 01

Taiwan
T: +886 2 2747 7278 F: +886 2 2747 7270

Turkey – Istanbul
T: +90 216 464 36 66 F: +90 216 464 36 67

UK (Freephone)
T: 0844 800 4512 F: 0844 800 4513

Uruguay
T: +59 82 900 6212 F: +59 82 900 6212

Venezuela
T: +58 212 959 6275 F: +58 212 959 5328

Boston Scientific

Delivering what's next.™

www.bostonscientific.com
www.bostonscientific-international.com

RCS Nanterre B420 668 420
© 2008 Boston Scientific Corporation
or its affiliates. All rights reserved.
DINEND2171EB

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for the use only in countries with applicable health authority product registrations.
PSST 4076 Printed in France by cwcreativeservices.